

Insulin Pen Devices

Which device matches which insulin?

With drug companies producing several types of insulins it is often difficult to work out which insulin cartridge fits with which insulin pen.

Non disposable insulin devices must be compatible with the insulin product they were designed for. Some insulin cartridges will simply not fit into the wrong device or if they do fit but it has not been tested in that device, an accurate dose delivery can not be guaranteed. Most companies produce the pens in different colours so those using more than one type of insulin can easily differentiate between the pen with short acting insulin from that with a longer acting insulin.

NOVO NORDISK

The below Novo Nordisk insulin cartridges can be used with these Novo Nordisk non disposable pens:

NOVO NORDISK Insulin Cartridges

- NovoRapid®
- Actrapid®
- Fiasp®
- NovoMix 30®
- Ryzodeg®
- Mixtard® 30/70 or 50/50
- Protophane®
- Levemir®

NOVO NORDISK Non Disposable Pens

- NovoPen 5
(can be dialed up to 60 units in 1 unit increments)




- NovoPen Echo
(can be dialed in ½ unit increments up to 30 units)
This pen has a memory is designed to assist people with remembering when insulin was last injected.



The **Innolet** is another **disposable** device from Novo Nordisk. The Innolet can be dialed in 1 unit increments up to 50 units.
Insulins available in the innolet include:

- Mixtard 30/70®
- Protophane®

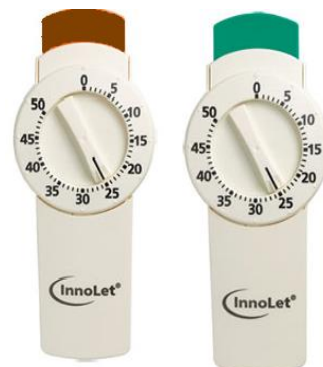

Darwin

Shop 1-3 Tiwi Place,
Tiwi NT 0810
PO Box 40113,
Casuarina NT 0811
Phone: 08 8927 8488
Fax: 08 8927 8515
E: info@healthylivingnt.org.au

Alice Springs

Jock Nelson Centre,
7/16 Hartley Street,
Alice Springs NT 0870
Phone: 08 8952 8000
Fax: 08 8952 7000
E: alicespings@
healthylivingnt.org.au

www.healthylivingnt.org.au
ABN 11 374 693 055

Healthy Living NT
is the trading name of the
Diabetes Association
of the Northern Territory
Incorporated.

Healthy Living NT
is the registered NT
licence holder for
Life. Be in it.



ELI LILY

The below Eli Lilly insulin cartridges can be used with these Eli Lilly non disposable pens:

ELI LILY Insulin Cartridges

- Humalog®
- Humalog Mix® 25/Humalog Mix® 50
- Humulin

ELI LILY Non Disposable Pens

- HumaPen® Luxura

(can be dialed up in 1 unit increments up to 60 units)

- HumaPen® Luxura HD

(can be dialed up in ½ unit increments up to 30 units)



SANOFI AVENTIS

The below Sanofi Aventis insulin cartridges can be used with these Sanofi Aventis non disposable pens:

SANOFI AVENTIS Cartridges

- Lantus®
- Apidra®
- Optisulin®



SANOFI AVENTIS Non Disposable Pens

- AllStar
- AllStar Pro
- JuniorStar
- ClikStar

(can be dialed up in 1 unit increments up to 80 units)



SANOFI AVENTIS Toujeo

- Toujeo® has a different concentration of insulin and only comes in pre-filled pens



MYLAN

- Semglee® – only comes in pre-filled pens



Note: All BD, PenFine and Novofine needles are compatible with all insulin pens and devices, disposable and non disposable.